


BIOCURVE BROAD

SCREW HOLE SEQUENCE

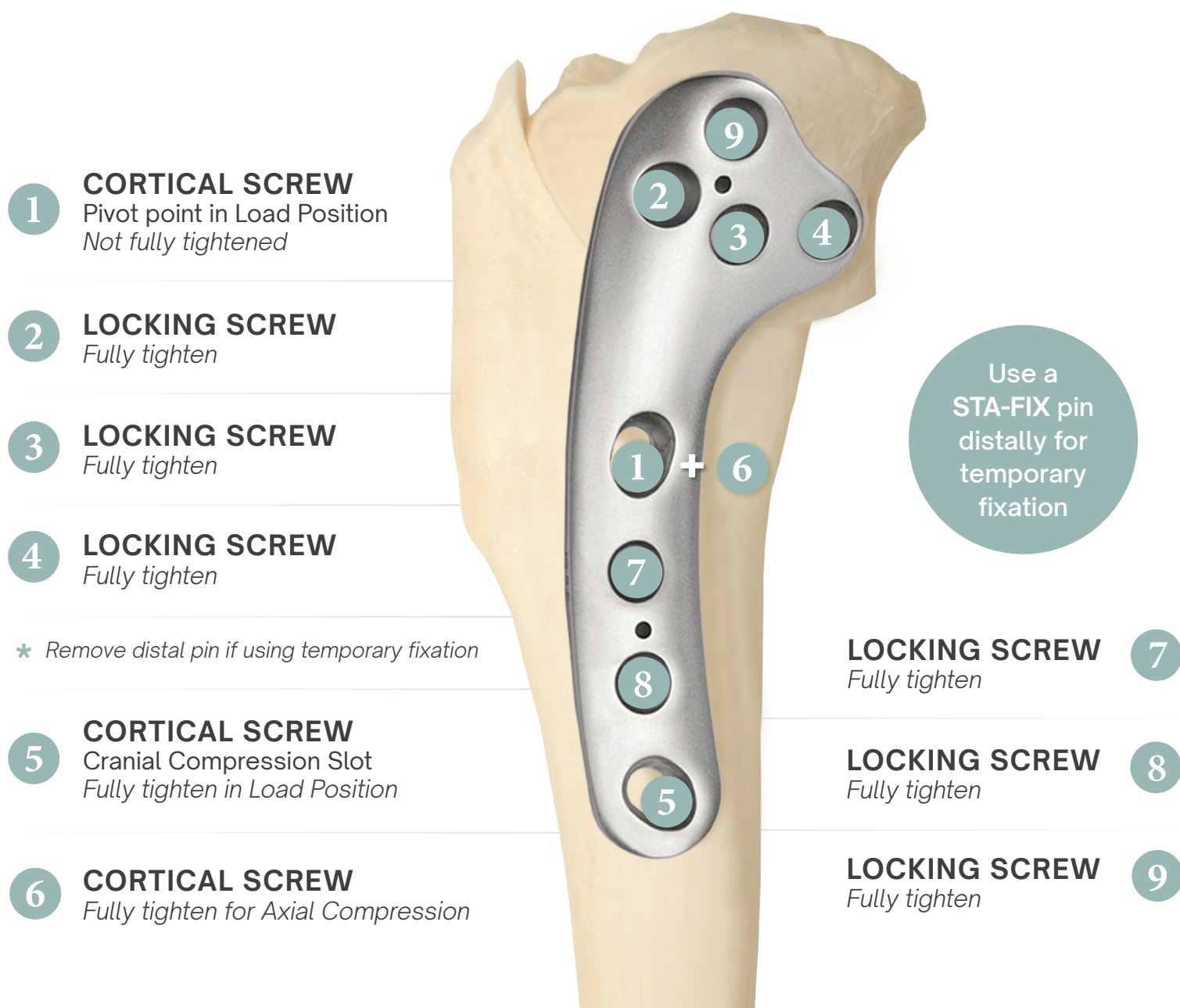

Note: The screw sequence is a recommended guide only.
In certain situations, this exact sequence may not be able to be executed.

©2024 Movora

BIOCURVE JUMBO

SCREW HOLE SEQUENCE

1 CORTICAL SCREW
Pivot point in Load Position
Not fully tightened


2 LOCKING SCREW
Fully tighten

3 LOCKING SCREW
Fully tighten

4 LOCKING SCREW
Fully tighten

* Remove distal pin if using temporary fixation

5 CORTICAL SCREW
Cranial Compression Slot
Fully tighten in Load Position

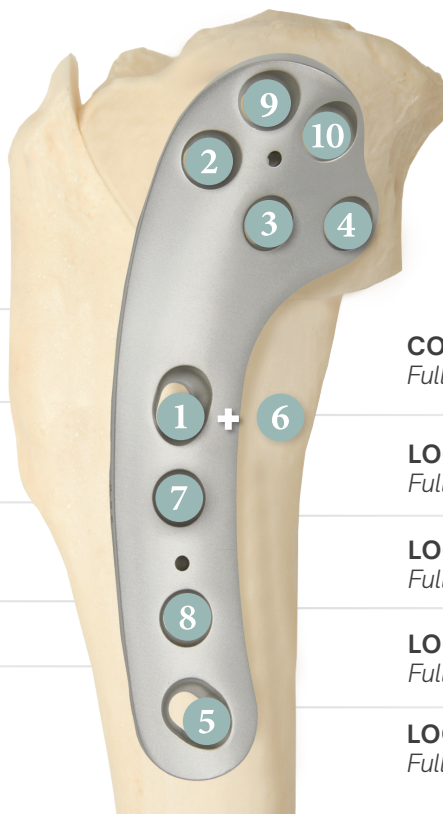

CORTICAL SCREW
Fully tighten for Axial Compression

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

Note: The screw sequence is a recommended guide only.
In certain situations, this exact sequence may not be able to be executed.

BIOCURVE JUMBO PLUS

SCREW HOLE SEQUENCE

1 CORTICAL SCREW
Pivot point in Load Position
Not fully tightened


2 LOCKING SCREW
Fully tighten

3 LOCKING SCREW
Fully tighten

4 LOCKING SCREW
Fully tighten

* Remove distal pin if using temporary fixation

5 CORTICAL SCREW
Cranial Compression Slot
Fully tighten in Load Position

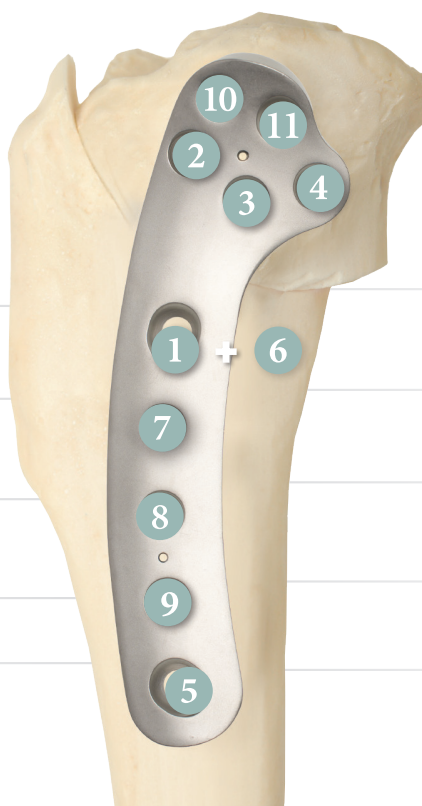

CORTICAL SCREW
Fully tighten for Axial Compression

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

LOCKING SCREW
Fully tighten

MMKT-SS-0056-A